

**Custom Laser
Engraving
Available!**


Flopac™

Fully integrated and self-contained Flopac offers Flotec's legendary durability, dependability, flexibility and safety in a unique all-in-one lightweight system incorporating several important design and safety features not found elsewhere.

- ✓ **Custom 0-3000psi and 0-4000psi External Pressure Contents Gauges**
- ✓ **Manual External Pressure Shut Off Knob**
- ✓ **Anti-Blowout Design with No Snap Rings**
- ✓ **24 Different Flow Ranges Available**
- ✓ **Full-Time Burst Disc Protection**
- ✓ **30 Micron Inlet Filtration**
- ✓ **All Metal Construction**
- ✓ **Outstanding Impact Resistance**
- ✓ **All Brass High Pressure Conduit**
- ✓ **Unique Halo' Style Handle Design (pictured)**
- ✓ **MRI Compatible and Safe**
- ✓ **No Metal-to-Metal Contact Points Along Internal Gas Flow Pathway**
- ✓ **None, One or Two 50psi DISS-1240 Check Valves**
- ✓ **Meets EN738-3 Specifications**
- ✓ **5-Year Warranty**

**0-3000 and 0-4000psi Pressure
Gauge Models Available**

**High Pressure Aluminum
and
Composite Cylinder Ready!**

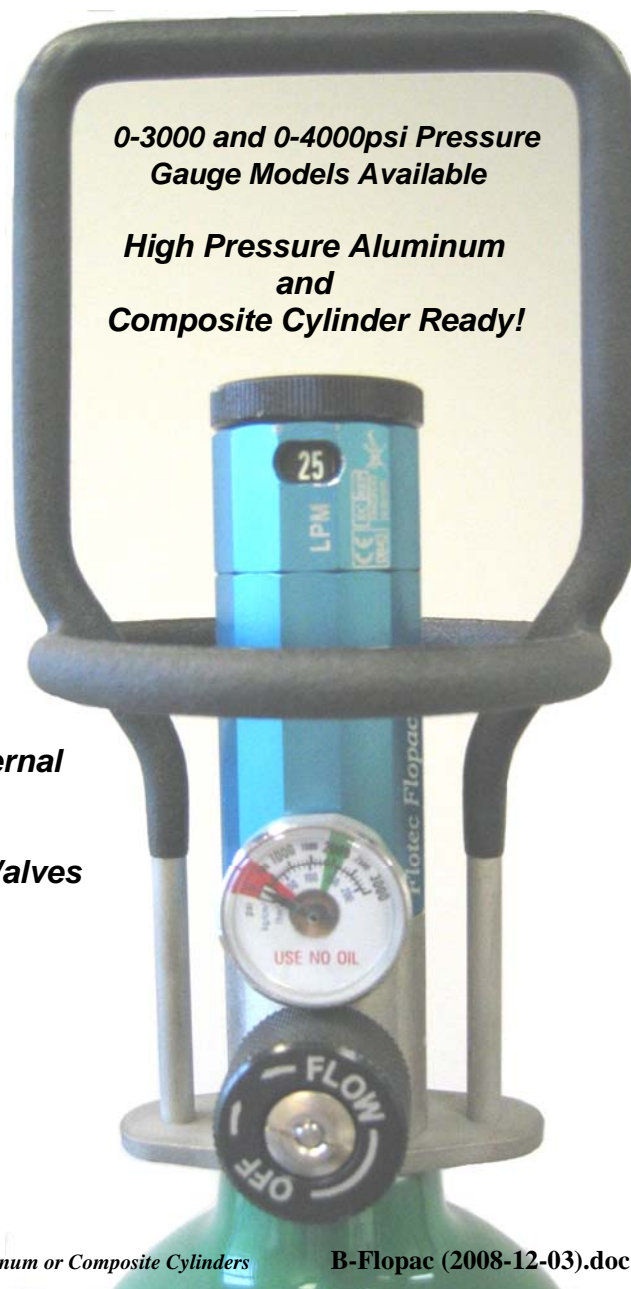



All Metal Construction (No Plastic Parts)

“Night Glow” Glow in the Dark Flow Numeral Decals

Solid Stainless-Steel ‘Halo’ Style Protective Handle

Distinctive Dodecagon (12 Sided) Body Style with 10 Different Anodized Colors Available

Custom Laser Text and Logo Engraving available

Interlocking Machined Threaded O-Ring Module Seals (Anti-Blowout Design using No Snap Rings)

Custom 0-3000psi or 0-4000psi External Contents Pressure Gauge

Safety Cylinder Pressure Shutoff Knob (EN738-3 Requirement)



Available 45° Angled Handle for Smaller Footprint EMS O₂ Responder Bags



Hose Barb or DISS-1240 Outlet Connection

30 Micron Stainless-Steel Sintered Inlet Filtration

Plain, One or Two 50psi DISS-1240 Check-Valve Outlets

Tri-Metal Plated Brass High-Pressure Connector Module with Oligodynamic Anti-Germicidal Characteristics

Standard CGA-540 Fill Port (Other Gases Available)

Aluminum Cylinder Sizes: C, D, Jumbo D, E

Composite Carbon Fiber Sizes: B, C, D equivalents

High Pressure Aluminum Sizes: B, C, D, E

ISO 9001:2000 & ISO 13485:2003 Certified